General Safety Information

INTENDED USE

Moisture Vue™ On-Demand Moisture Chamber is an eye shield that attaches to ophthalmic frames and functions as a moisture chamber to help maintain eye moisture for young adults and adults who wear eyeglasses.

DIRECTIONS FOR USE

(See Figure on the Right)

- 1) Clean glasses surface with appropriate lens cleaner.
- 2) Remove protective film on device to expose adhesive.
- 3) Affix adhesive side of device to outer edge of glasses frame or outer edge of lenses. Once device is in appropriate position, firmly pressed the device to the glasses.
- 4) Check for gaps. For maximum moisture retention, ensure the seal sits flush against the skin without gaps. Reposition device to your comfort by removing the device from the glasses and reapplying the device if necessary.
- 5) Repeat process on other side of glasses.

STORAGE INSTRUCTIONS

- Place device in the provided storage case in between uses.
- Store chambers in room-temperature, low humidity conditions.
- Avoid leaving chambers in direct sunlight or extreme heat (e.g., car dashboards).
- Keep away from sharp objects that may damage seals.

CLEANING & REUSABILITY INSTRUCTIONS

- Gently wipe the chambers with a dry cloth between uses.
- Avoid getting water, soap, abrasive cleaners, solvents, or alcohol-based solutions on the adhesive surface to maintain longevity.
- This device is a limited use device.
- Replace device once adhesive no longer provides a secure hold and dispose of device in accordance with local recycling regulations.
- This device must be stored between uses as indicated. Inspect for damage and adhesive hold prior to re-use.

CONTRAINDICATIONS

- Non-Sterile: Do not use on broken/irritated skin or if you have an untreated eye infection.
- Not intended for use as industrial safety eyewear or impact protection.
- Do not use if you have known allergies to silicone, adhesives, or frame materials.

WARNINGS

- Do not use while driving, cycling, operating heavy machinery, or performing tasks that require full peripheral vision. This device may restrict peripheral vision or fog unexpectedly.
- Discontinue use immediately if you experience pain, dizziness, allergic reaction, or visual distortion.
- Do not use in environments requiring certified eye protection (e.g., industrial safety, impact hazards) — these chambers are not safety-rated and are not a replacement for certified safety glasses.
- Small part choking risk, keep out of reach of children under 5yr.

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
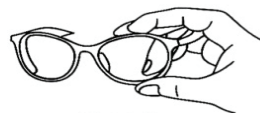
Step 5
Repeat on the other side

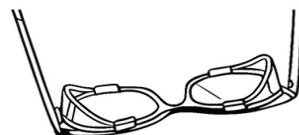CAUTIONS

- Avoid exposure to extreme heat (e.g., car dashboards, ovens) — silicone and adhesives may warp or degrade.
 - Do not clean with harsh chemicals (acetone, bleach, alcohol) — may damage lenses, coatings, or seals.
 - Do not submerge in water.
- Handle gently when inserting/removing chambers — excessive force may damage the frame, lens, or chamber
- Handle the device with care when inserting or removing the compartments: excessive force could damage the frame, lenses, or compartments.
 - Use the device at your own risk. Briter Eyez accepts no liability for damage caused to the glasses.

PRECAUTIONS

- Ensure a proper fit: the eyewear should be comfortable and must not exert excessive pressure on the face.
- Replace the device regularly to prevent oils or debris from compromising the seal.
- Check the integrity of the adhesive/seal before each use; replace it if it is damaged or peeling.
- Use only as directed: this device is intended to relieve symptoms of dry eyes and maintain eye moisture, not to treat other medical conditions.
- Consult your ophthalmologist before using this product* if you have an eye disease, have undergone recent surgery, or suffer from a chronic condition.
- Remove the eyewear periodically to rest your eyes and prevent skin irritation caused by prolonged contact.

For current patent and trademark information, please visit our website at www.moisturevue.com

For questions, feedback, or to report a complaint, please contact us at support@moisturevue.com

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